

**KAMIYA BIOMEDICAL COMPANY**

# Dog IgM ELISA

**For the quantitative determination of IgM  
in dog biological samples.**

**Cat. No. KT-458**

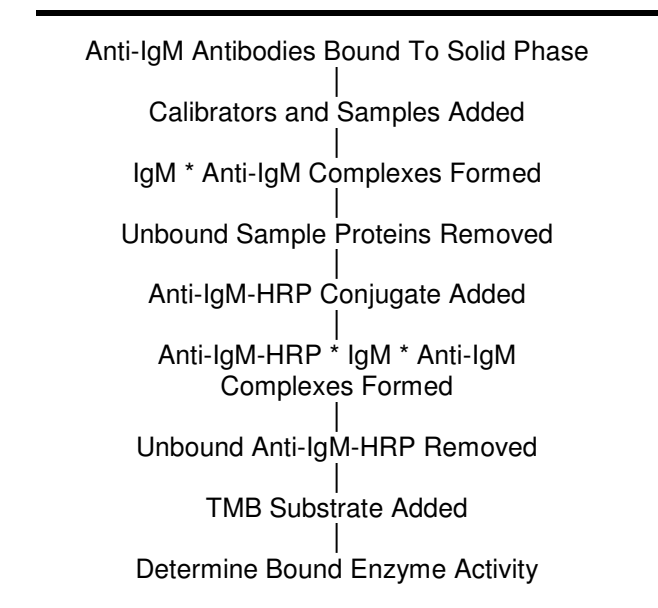
**For Research Use Only.**

**PRODUCT INFORMATION****Dog IgM ELISA**  
**Cat. No. KT-458****INTENDED USE**

The Dog IgM ELISA is a highly sensitive two-site enzyme-linked immunoassay (ELISA) for the quantitative determination of IgM in dog biological samples. For research use only.

**PRINCIPLE**

The principle of the double antibody sandwich ELISA is represented in Figure 1. In this assay the IgM present in sample reacts with the anti-IgM antibodies which have been adsorbed to the surface of polystyrene microtitre wells. After the removal of unbound proteins by washing, anti-IgM antibodies conjugated with horseradish peroxidase (HRP), are added. These enzyme-labeled antibodies form complexes with the previously bound IgM. Following another washing step, the enzyme bound to the immunosorbent is assayed by the addition of a chromogenic substrate, 3,3',5,5'-tetramethylbenzidine (TMB). The quantity of bound enzyme varies directly with the concentration of IgM in the sample tested; thus, the absorbance, at 450 nm, is a measure of the concentration of IgM in the test sample. The quantity of IgM in the sample can be interpolated from the calibration curve constructed from the calibrators, and corrected for sample dilution.

**Figure 1.****COMPONENTS**

1. Diluent Concentrate  
One bottle containing 50 mL of a 5X concentrated diluent running buffer.
2. Wash Solution Concentrate  
One bottle containing 50 mL of a 20X concentrated wash solution.
3. Enzyme-Antibody Conjugate Concentrate  
One vial containing 150 µL of a 100X concentrated affinity-purified anti-dog IgM antibody conjugated with HRP in a stabilizing buffer.
4. TMB Substrate Solution  
One vial containing 12 mL of TMB and hydrogen peroxide in citric acid buffer at pH 3.3.
5. Stop Solution  
One vial containing 12 mL of 0.3 M sulfuric acid. WARNING: Avoid contact with skin.

6. Microtiter Plate  
Twelve removable eight-well strips in well holder frame. Wells are coated with affinity-purified anti-dog IgM.
7. Dog IgM Calibrator  
One vial containing 0.1 mL of a Dog IgM Calibrator.
8. Positive Control  
One vial containing 50  $\mu$ L of serum with 0.1% sodium azide. See the Control Certificate for the concentration.

## MATERIALS REQUIRED BUT NOT PROVIDED

- Precision pipettes (2  $\mu$ L to 200  $\mu$ L) for making and dispensing dilutions.
- Test tubes
- Microplate washer/aspirator
- Distilled or de-ionized H<sub>2</sub>O
- Microplate reader
- Assorted glassware for the preparation of reagents and buffer solutions
- Timer
- Vortex mixer

## PRECAUTIONS

1. Read the instructions carefully before beginning the assay.
2. This kit is for research use only.
3. Great care has been taken to ensure the quality and reliability of this product. However, it is possible that in certain cases, unusual results may be obtained due to high levels of interfering factors.
4. Preservatives
  - Positive Control contains 0.1% sodium azide.
5. No additives or preservatives are necessary to maintain the integrity of the specimen. Avoid azide contamination.
6. Azide and thimerosal at concentrations higher than 0.1% inhibit the enzyme reaction.
7. Other precautions:
  - Do not interchange kit components from different lots.
  - Do not use kit components beyond the expiration date.
  - Protect reagents from direct sunlight.
  - Do not pipette by mouth.
  - Do not eat, drink, smoke or apply cosmetics where reagents are used.
  - Avoid all contact with the reagents by using gloves.
  - Stop solution contains diluted sulfuric acid. Irritation to eyes and skin is possible. Flush with water after contact.

## REAGENT PREPARATION

1. Diluent Concentrate  
The Diluent Solution supplied is a 5X Concentrate and must be diluted 1:5 with distilled or deionized water.
2. Wash Solution Concentrate  
The Wash Solution supplied is a 20X concentrate and must be diluted 1:20 with distilled or de-ionized water. Crystal formation in the concentrate is not uncommon when storage temperatures are low. Warming of the concentrate to 30-35°C before dilution can dissolve crystals.
3. Enzyme-Antibody Conjugate Concentrate  
Calculate the required amount of working conjugate solution for each microtiter plate test strip by adding 10  $\mu$ L Enzyme-Antibody Conjugate to 990  $\mu$ L of 1X Diluent for each test strip to be used for testing. Mix uniformly, but gently. Avoid foaming.
4. TMB Substrate Solution  
Ready to use as supplied.
5. Stop Solution  
Ready to use as supplied.
6. Microtiter Plate  
Ready to use as supplied. Unseal Microtiter Pouch and remove plate from pouch. Remove all strips and wells that will not be used in the assay and place back in pouch and re-seal along with desiccant.

## 7. Dog IgM Calibrator

The calibrator is now at a concentration of 1.7 mg/mL. Dog IgM calibrators need to be prepared immediately prior to use (see chart below). Mix well between each step. Avoid foaming.

Calibrator	Concentration (ng/mL)	Calibrator Volume added to 1X Diluent	Volume of 1X Diluent
7	17,000	10 $\mu$ L Dog IgM Calibrator	990 $\mu$ L
6	800	40 $\mu$ L Calibrator 0	810 $\mu$ L
5	400	250 $\mu$ L Calibrator 1	250 $\mu$ L
4	200	250 $\mu$ L Calibrator 2	250 $\mu$ L
3	100	250 $\mu$ L Calibrator 3	250 $\mu$ L
2	50	250 $\mu$ L Calibrator 4	250 $\mu$ L
1	25	250 $\mu$ L Calibrator 5	250 $\mu$ L
0	0	0	500 $\mu$ L

## 8. Positive Control

The concentration and recommended dilution are provided on the Control Certificate. Before use, briefly centrifuge the Positive Control to allow all of the liquid to collect in the bottom of the vial.

**STORAGE AND STABILITY**

## 1. Complete Kit

The expiration date for the kit is stated on the outer label. The recommended storage temperature is 4°C. **Note: See long-term storage recommendations below for the Positive Control.**

## 2. Diluent

The 5X Diluent Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from date of preparation. Both solutions should be stored at 4°C.

## 3. Wash Solution

The 20X Wash Solution Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions can be stored at room temperature (RT, 16-25°C) or at 4°C.

## 4. Enzyme-Antibody Conjugate

100X concentrated horseradish peroxidase anti-dog IgM conjugate should be stored at 4°C and diluted immediately prior to use. The working conjugate solution is stable for up to eight hours.

## 5. TMB Substrate Solution

The TMB Substrate Solution should be stored at 4°C and is stable until the expiration date.

## 6. Stop Solution

The Stop Solution should be stored at 4°C and is stable until the expiration date.

## 7. Microtiter Plate

Anti-dog IgM coated wells are stable until the expiration date and should be stored at 4°C in the sealed foil pouch with a desiccant pack.

## 8. Dog IgM Calibrator

The Dog IgM Calibrator should be stored at 4°C. The working calibrator solutions should be prepared immediately prior to use and are stable for up to eight hours.

## 9. Positive Control

For storage longer than 7 days keep frozen until the expiration date. Storage less than 7 days can be at 4°C. Avoid multiple freeze/thaw cycles.

**INDICATIONS OF INSTABILITY**

If the test is performing correctly, the results observed with the calibrator solutions should be within 20% of the expected values.

## SPECIMEN COLLECTION AND HANDLING

Blood should be collected by venipuncture and the serum separated from the cells, after clot formation, by centrifugation. For plasma samples, blood should be collected into a container with an anticoagulant and then centrifuged. Care should be taken to minimize hemolysis, excessive hemolysis can impact your results. Assay immediately or aliquot and store samples at -20°C. Avoid repeated freezing/thawing.

For any sample that might contain pathogens, care must be taken to prevent contact with open wounds. No additives or preservatives are necessary to maintain the integrity of the specimen. Avoid azide contamination.

## ASSAY PROTOCOL

### Dilution of Serum Samples

Due to the high-sensitive nature of the assay, each test sample should be diluted before use for a normal assay. A 1:10,000 dilution is appropriate for most serum samples. For absolute quantification of samples that yield results outside the range of the calibration curve, a lesser or greater dilution might be required.

To prepare a 1:10,000 dilution of sample, transfer 5  $\mu$ L of sample to 495  $\mu$ L of 1X Diluent. This gives you a 1:100 dilution. Next, dilute the 1:100 samples by transferring 5  $\mu$ L to 495  $\mu$ L of 1X Diluent. You now have a 1:10,000 dilution of your sample. Mix thoroughly at each step.

### Procedure

1. Bring all reagents to RT before use.
2. Pipette 100  $\mu$ L of
  - Calibrator 1 (800 ng/mL) in duplicate
  - Calibrator 2 (400 ng/mL) in duplicate
  - Calibrator 3 (200 ng/mL) in duplicate
  - Calibrator 4 (100 ng/mL) in duplicate
  - Calibrator 5 (50 ng/mL) in duplicate
  - Calibrator 6 (25 ng/mL) in duplicate
  - Calibrator 0 (0 ng/mL) in duplicate
3. Pipette 100  $\mu$ L of diluted Positive Control (in duplicate) into pre designated wells.
4. Pipette 100  $\mu$ L of diluted sample (in duplicate) into pre designated wells.
5. Incubate the Microtiter Plate at 22°C (RT) for thirty (30  $\pm$  2) minutes. Keep plate covered and level during incubation.
6. Following incubation, aspirate the contents of the wells.
7. Completely fill each well with appropriately diluted Wash Solution and aspirate. Repeat three times, for a total of four washes. If washing manually: completely fill wells with wash buffer, invert the plate then pour/shake out the contents in a waste container. Follow this by sharply striking the wells on absorbent paper to remove residual buffer. Repeat three times for a total of four washes.
8. Pipette 100  $\mu$ L of appropriately diluted Enzyme-Antibody Conjugate to each well. Incubate at 22°C (RT) for thirty (30  $\pm$  2) minutes. Keep plate covered in the dark and level during incubation.
9. Wash and blot the wells as described in Steps 6 and 7.
10. Pipette 100  $\mu$ L of TMB Substrate Solution into each well.
11. Incubate in the dark at RT for precisely ten (10) minutes.
12. After ten minutes, add 100  $\mu$ L of Stop Solution to each well.
13. Determine the absorbance at 450 nm of the contents of each well. Zero the plate reader to air.

The absorbance of the final reaction mixture can be measured up to 2 hours after the addition of the Stop Solution. However, good laboratory practice dictates that the measurement be made as soon as possible.

## RESULTS

1. Subtract the average background value from the test values for each sample.
2. Using the results observed for the calibrators construct a calibration curve. The appropriate curve fit is that of a four-parameter logistics curve. A second order polynomial (quadratic) or other curve fits may also be used.
3. Interpolate test sample values from the calibration curve. Correct for sample dilution factor to arrive at IgM concentration in original sample.

## QUALITY CONTROL

In accord with good laboratory practice, the assays for total IgM require meticulous quality control. Each laboratory should use routine quality control procedures to establish inter- and intra-assay precision and performance characteristics.

## LIMITATION OF THE PROCEDURE

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the information contained in the package insert instructions and with adherence to good laboratory practice.
2. Factors that might affect the performance of the assay include proper instrument function, cleanliness of glassware, quality of distilled or de-ionized water, and accuracy of reagent and sample pipettings, washing technique, incubation time or temperature.

### **FOR RESEARCH USE ONLY**

#### **KAMIYA BIOMEDICAL COMPANY**

12779 Gateway Drive, Seattle, WA 98168  
Tel: (206) 575-8068 Fax: (206) 575-8094  
Email: LifeScience@k-assay.com  
www.k-assay.com